THOMSON

Hydropulseur TH3HYDRO

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THOMSON

TH3HYDRO

WATERPIK
HYDROPULSEUR
WATERFLOSSER
IRRIGADOR
IRRIGADOR

 OPERATING INSTRUCTIONS

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 HANDLEIDING

 MANUAL DE INSTRUCCIONES

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TABLE OF CONTENTS

IMPORTANT SAFETY INSTRUCTIONS.....	3
PRODUCT OVERVIEW.....	6
CHARGING THE APPLIANCE.....	6
USING THE APPLIANCE.....	7
CLEANING AND MAINTENANCE.....	8
SPECIFICATIONS.....	8
DISPOSAL.....	9

IMPORTANT SAFETY INSTRUCTIONS

PLEASE READ THE FOLLOWING SAFETY WARNINGS AND THE INSTRUCTIONS CAREFULLY BEFORE USING THE APPLIANCE AND KEEP FOR FUTURE REFERENCE.

Warnings

- Read the instructions.
- This appliance is intended for domestic household use only and should not be used for any other purpose or in any other application, such as for non-domestic use or in a commercial environment.
- This oral irrigator can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- This oral irrigator shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- The supply terminals are not to be short-circuited.
- Attention should be drawn to the environmental aspects of battery disposal. Don't throw used batteries in dustbin. Please contact your retailer in order to protect the environment.
- The batteries (batteries installed) shall not be exposed to excessive heat such as sunshine, fire or the like.

- If battery leakage occurs, avoid contact with skin. Isolate the leaking battery in a sealed plastic bag and put it in the scrap following the instructions environmental protection. If contact with the skin, mucous membranes or eyes, rinse thoroughly with water and contact your doctor or ophthalmologist.
- Regarding the detail on how to remove the battery, refer to the section “Remove the Battery” in pages EN-5 & EN-6.
-  Detachable supply unit
- Use only with BI12TA-050060-Bdv supply unit.
- This appliance contains batteries that are non-replaceable.
- **WARNING:** For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. (Adaptor model number: BI12TA-050060-Bdv, Made by: Blue Iron Holdings Limited)
- The external flexible cable of the transformer cannot be replaced, if the cord of this transformer is damaged, it is necessary to dispose the transformer and replaced by a transformer of the same model.

Safety Information

- Toddlers, infant, people who are unable to operate and people with weak oral senses should not use this appliance.
- Do not use for any purpose other than oral cleaning.
- Do not fill the tank with water warmer than 40°C.
- Do not use the appliance while bathing.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not reach for an appliance that has fallen into water. Unplug immediately.


- Keep the cord away from heated surfaces.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can become lodged in the throat.
- Certain types of mouthwash may damage the appliance by causing the housing or water tank to crack. Therefore, do not use mouthwash in the water tank.
- People who experience tooth or gum pain may not be able to use this appliance and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums.
- Never alter the appliance. Also, do not try to disassemble or repair it.

Remove the Battery

The appliance is supplied with one rechargeable Lithium 3.7V 1400mAh battery incorporated in the appliance which is not accessible and cannot be replaced.

Therefore, when the appliance reaches the end of its useful life it must be disposed of at a recycling centre.


When disposing of the appliance, first remove the rechargeable battery from the housing by following the instructions below.

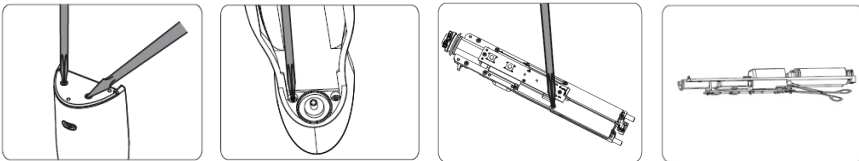1. Remove the 2 covers from the bottom of the housing and then remove the bottom.
2. Slide the water tank off the housing and unscrew the screws.
3. Put the inner part out of the appliance, then unscrew all screws fixed on the inner part.

4. Cut off the line and remove the battery from the inner part.

PRODUCT OVERVIEW

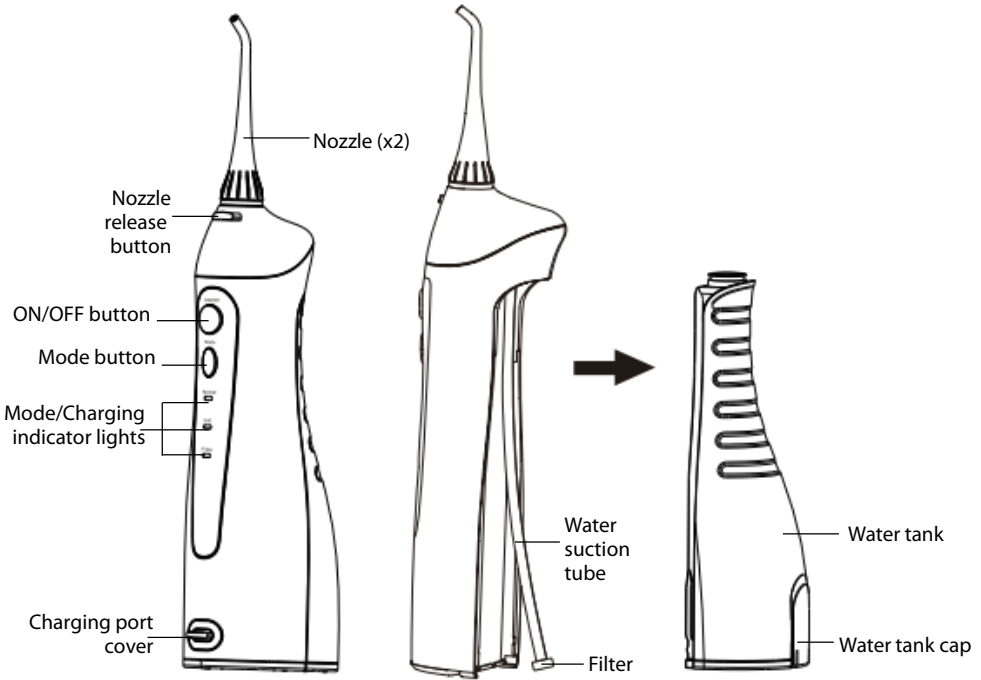
Note:

There may be some liquid remaining inside the main appliance (tank, water suction tube, etc.), but this is remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

CHARGING THE APPLIANCE

It takes about 4 hours to fully charge the appliance.
A full charge lasts for 30 to 50 minutes of use.

1. Remove the charging port cover.
2. Plug the supplied power adapter into the charging port.
3. Plug the adapter into the mains socket.
4. The CH indicator light will flash during charging and will go out when the appliance is fully charged.
5. Attach the charging port cover back to the appliance.

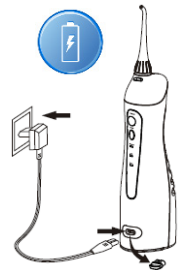
Note:

If the CH indicator light flashes during or after an irrigating, recharge the battery.

USING THE APPLIANCE

1. Attach the nozzle to the appliance.
2. Open the water tank cap, fill the tank with water which temperature should not exceed 40°C, then close the cap.

Or, remove the water tank from the appliance, hold the tank upright, fill the tank with water (temperature not exceed 40°C), then attach the tank onto the appliance.

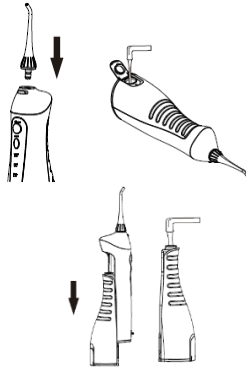3. Repeatedly press the Mode button to select your desired mode.
 - The corresponding mode indicator light will illuminate.

3 Different Modes


NORMAL

For normal cleaning.

SOFT


Massage & stimulate gums for enhanced circulation.

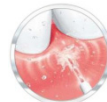PULSE

Clean deep between teeth and below the gum line.


4. Hold the appliance in an upright position with nozzle tip pointed toward your teeth.
5. Keep your mouth opened slightly so that the water can flow out.
6. Press the ON/OFF button to start irrigating.

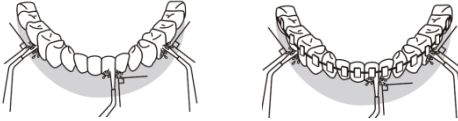- Direct the water jet flow at a 90-degree angle to the teeth and gums.
 - Direct the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue.
 - Clean between the teeth and around any braces, crowns or bridgework thoroughly.
7. Once irrigation is complete, press the ON/OFF button to turn off the appliance.


Note:

If you need to replace the nozzle, press the ON/OFF button to turn off the appliance and then press the nozzle release button to remove the nozzle. Then attach another nozzle to the appliance.

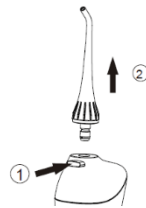After using the appliance

1. Open the water tank cap and dispose of the remaining water in the tank.
2. Press the ON/OFF button to turn on the appliance, and drain the remaining water from inside the appliance.
3. Press the ON/OFF button again to turn off the appliance, and close the water tank cap.
4. Wipe the appliance dry with a soft cloth.

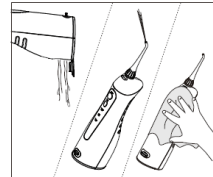To prevent the growth of bacteria, make sure there is no remaining water inside the appliance. If you are not going to use the appliance for some time, wipe and dry the appliance thoroughly before storing it.

Two-Minute Auto Timer

The appliance will shut off automatically after 2 minutes of use.
Two minutes auto-timer can avoid it still working after you finish cleaning.

CLEANING AND MAINTENANCE

IMPORTANT:

Turn off the appliance before cleaning.

Use only water or neutral cleaning agent to clean the appliance.

Do not wash using hot water over 50°C. Keep the appliance away from high temperature area and avoid direct sunlight.

Main Unit

Clean the main unit with a damp cloth.

Do not immerse the main unit in water.

Wipe the unit dry with a soft cloth.

Water Tank Cap

Wipe off any grime and dirt with a cloth.

Tank

Wash with water and wipe dry with a soft cloth.

Empty the tank when you are not going to use it for a long period.

Nozzle

Wash with water and wipe with a soft cloth.

Water Suction Tube

Wash with water and wipe with a soft cloth.

Do not bend, pull or twist the water suction tube.

SPECIFICATIONS

Adapter input: 100-240V~, 50/60Hz, 0.5A
output: 5V--- 600mA

Smart oral irrigator input: 5V--- 600mA

Thomson and the Thomson logo are trademarks of Technicolor (S.A.) or its affiliates and are used under license by Fnac Darty (S.A.).

DISPOSAL

As a responsible retailer we care about the environment. As such we urge you to follow the correct disposal procedure for the appliance and packaging materials. This will help conserve natural resources and ensure that it is recycled in a manner that protects health and the environment.

You must dispose of this appliance and its packaging according to local laws and regulations. Because this appliance contains electronic components, the appliance and its accessories must be disposed of separately from household waste when the appliance reaches its end of life.

Contact your local authority to learn about disposal and recycling. The appliance should be taken to your local collection point for recycling. Some collection points accept appliance free of charge.

We apologise for any inconvenience caused by minor inconsistencies in these instructions, which may occur as a result of product improvement and development.

Etablissements Darty & fils ©, 129 Avenue Gallieni, 93140 Bondy, France 13/10/2020

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/Modelnummer / Numero de modelo /

Numero do modelo:

We declare that the above product has been evaluated and deemed to comply with the requirements of the listed European Directives and standards:

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- 2014/30/EU
- 2014/35/EU
- 2009/125/EC
- 2011/65/EU & (EU) 2015/863

Title

EMC Directive (EMC)
Low Voltage Directive (LVD)
Eco design requirements for energy-related products (ErP)
RoHS

EMC:

EN 55014-1:2017
EN 55014-2:2015
EN 61000-3-2:2014
EN 61000-3-3:2013

LVD:

EN 60335-1:2012+A11:2014+A13:2017
EN 60335-2-52:2003+A1:2008+A11:2010
EN 62233:2008

ErP:

(EC) No 278/2009
(EU) No 2016/2282
EN 50563:2011+A1:2013

The person responsible for this declaration is:

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Name - Nom-: Hervé HEDOUIS

Position - Poste: General Manager Direct Sourcing

Signature - Signature:

A handwritten signature in blue ink, appearing to read 'Hedouis'.

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Bondy, July 15, 2020